AquaFloss Compact C300

EN - USER MANUAL

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HR - UPUTE ZA UPOTREBU

EN - USER MANUAL


BOX CONTENTS

- TrueLife AquaFloss Compact C300
- 3x Standard nozzle heads, 1x Dental Plaque nozzle head, 1x Orthodontic nozzle head, 1x Periodontal nozzle head
- USB charging cable
- Manual and safety instructions

PRODUCT DESCRIPTION

1. Nozzle head
2. Button to release nozzle head
3. On/Off button
4. Display
5. Mode change button
6. Handle
7. Water tank
8. Water tank cap
9. Charging Cable

Warning: *There may be residual water in the device. This is distilled water used during quality control testing.*

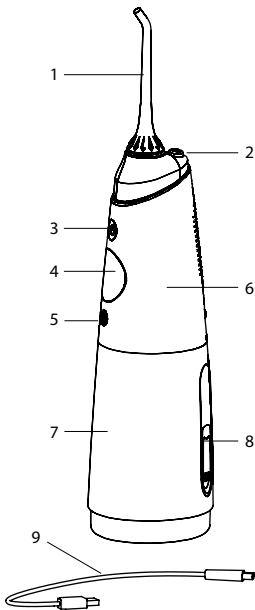INSTRUCTIONS FOR USE

Before the first use, remove all packaging material.

If the water flosser has been stored in a low temperature environment during transport, leave it at room temperature for at least two hours before use.

1. Attach the head to the water flosser. You will hear a clicking sound when it is correctly attached.
2. Remove the water tank cap and fill the tank with water; the temperature of the water should not exceed 40 °C. The capacity of the water tank is 300 ml. Then replace the cap, ensuring that water does not spill.
3. Press the „mode change“ button and select the desired flossing mode. The mode you have just selected will be shown on the display. You can also change the mode while flossing. You can choose between 4 modes:

Soft mode – suitable for people with sensitive teeth and for people starting out with oral shower. Pressure range: 20–51 PSI

Medium mode – suitable for everyday brushing. Pressure range: 50–80 PSI

Strong mode – suitable for deep cleaning. Pressure range: 70–110 PSI

Pulse mode – suitable for massaging teeth and gums. Pressure range: 70–100 PSI

DIY mode – adjust the pressure to suit your needs. Long press the mode change button to enter DIY mode. A scale will appear on the display and gradually go from minimum pressure to maximum pressure. Press the mode change button again to select the water pressure as required. Pressure range: 8–110 PS

4. Leave your mouth slightly open to allow water to rinse from your mouth into the sink. Press the On / Off button to start the flosser. Direct the water jet at a 90° angle to the teeth and gums. Slowly run the stream of water around your teeth. Direct the stream of water onto the gums and into the space between the teeth. If you have braces, you can direct the water between the joints of the braces.
5. When you have finished flossing, switch off the device using the On / Off button. Open the water tank cap and pour out any remaining water. Switch the device on again with the On / Off button to drain the water from the pump. Press the On / Off button to turn off the device and close the water tank cap.

Warning: To prevent the build-up of bacteria, ensure that no water remains in the device. If you do not intend to use the device for a longer time, make sure that it is thoroughly dried before storing it.

CHARGING

The water flosser is not supplied with a fully charged battery. Before using for the first time, fully charge the water flosser.

Plug one end of the USB cable into the charging socket at the rear of the water flosser. Plug the other end of the USB cable into an adapter (adapter not included). Insert the adapter into the plug socket.

It takes approximately 4 hours to fully charge the water flosser.

Note: To ensure maximum battery life, we recommend charging the water flosser once a week.

INSERTING AND REPLACE THE NOZZLE HEAD

Insert the nozzle head into the hole at the top of the water flosser. You will hear a click when the nozzle head has been correctly attached.

To remove the nozzle head from the main body, press and hold the release button and then detach the nozzle head from the water flosser.

The nozzle head rotates 360°.

Warning: Only insert and remove the nozzle head when the water flosser is switched off.

Note: Nozzle heads are items of personal hygiene and are not covered by warranty.

MAINTENANCE

Use only water and neutral detergents to clean the device. Never use caustic or abrasive cleaners as they may damage the device. Do not use water whose temperature exceeds 50 °C. Keep the device away from high temperatures and avoid direct sunlight.

1. The device

Remove dirt by wiping with a towel or cloth with liquid soap and water, then ensure all water is rinsed off and wipe dry.

Do not immerse the device in water.

Warning: Water must not enter the USB socket.

2. Water tank

Wash the tank with water and wipe with a dry cloth.

Dry the tank thoroughly if you do not intend to use it for a longer time.

Note: For greater ease of cleaning, the water tank can be removed from the body of the water flosser.

3. Nozzle head

Wash the nozzle head with water and wipe with a cloth.

Replace the nozzle head regularly, at least once every three months, or more often depending on wear.

TROUBLESHOOTING

Problem	Possible cause	Solution
Water pressure is very low	The nozzle head is unclean / battery is low	Replace the nozzle head with a new one / recharge the water flosser
No water at all emits from the water flosser	Water tank is empty / The device is tilted at the wrong angle	Refill the water tank / Move the device to the correct position (i.e., vertical)
Nozzle head cracks	Nozzle head has been used for too long	Replace the nozzle head with a new one
Water tank leaks	Water tank cap is damaged	Contact the manufacturer to replace the water tank

SPECIFICATIONS

Energy consumption	5 W
Voltage	5 V
Battery capacity	1400 mAh
Water resistance	IPx7
Pressure range	8–110 PSI
Reservoir volume	300 ml
Weight	239 g
Dimensions	63×68×304 mm

PRODUCT LIABILITY

TrueLife products are guaranteed to be free from manufacturing defects. In the event of such defects, elem6 s.r.o. will repair or replace the defective part or product with a comparable part or product at its sole discretion. The warranty will be deemed void if the device is altered or tampered with. The warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss, or other claims related to or resulting from the use of this product. elem6 s.r.o. shall in no event be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Return or exchange of spare parts and / or products may be subject to shipping, handling, and replacement or replenishment fees.

This product should be stored out of the reach of children. Use by children should take place only under the direct supervision of an adult.

SAFETY GUIDELINES WHEN OPERATING THE DEVICE

1. Read the enclosed safety instructions carefully.
2. The latest version of the manual, including instructions for operating the device, can be downloaded at: www.truelife.eu.
3. Printing errors and changes of information reserved.

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