FitScale W3

EN – INSTRUCTIONS FOR USE


The Truelife FitScale W3 smart wireless scales measures body weight using bioelectric impedance analysis (BIA). BIA technology is used to measure body fat and gather data on the body's overall composition. The Truelife FitScale W3 measures weight, fat, water, bone mass, muscle mass, calories, and body mass index (BMI). The scale is for home use only.

To achieve consistent and accurate readings, please use at the same time the day, under approximately the same conditions, if possible. Avoid weighing in the following situations:

1. During or immediately after a hard workout
2. Just after a shower, bath or sauna
3. After high alcohol consumption
4. During or immediately after a meal
5. During a fever

The unit of weight measurement can be changed from kilograms (kg) to pounds (lb) by pressing the button at the bottom of the scale. Height units (metre / foot) will also be changed.


Position the scales correctly on a flat surface or a solid surface. The weight should show 0.0 kg before use.

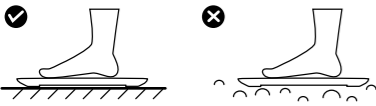Profile of scales on level ground and on uneven ground

To measure body fat and other body parameters correctly: stand upright with even weight distribution on both feet. The scale automatically switches on when you stand on it. Weigh yourself barefoot and stand straight. Weighing yourself whilst wearing socks will give inaccurate results. Dirty feet can also affect the readings.

After measuring your body weight, your scale will measure body parameters using BIA which will appear on the display. You can check the results within a few seconds, and the results are repeated on the display.

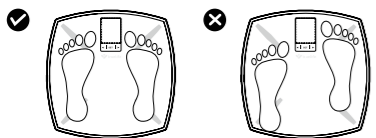Graphical representation of the feet on the weighing area

Settings:

Truelife FitScale W3 personal scale supports multiple users (P1-P8). Follow the instructions below to set up your personal profile.

1. Profile Settings

Hold your finger on SET to enter Settings. The system asks you to select the profile number (P1-P8). The icon „P1“ blinks as shown below. To select a different profile, press the keys ◀ or ▶. Press the „SET“ button to confirm the setting.

2. Gender Settings

After selecting the profile number, the system automatically moves on to setting the gender. Like the image below the icon ♂ flashes. The user can hold the button ◀ or ▶ to select gender. Press the „SET“ button to confirm the setting.

Athlete Mode: Athlete is defined as an adult who is engaged in intense physical activity about 12 hours a week and has a heart rate at rest of about 60 beats per minute or less.

3. Height Settings


After you confirm your gender, the system will guide you to the height selection. As shown in the figure below, the number flashes on the display. By holding down the button ◀ or ▶ the user sets their height. (Range 100 cm to 220 cm). Press the „SET“ button to confirm the setting.

4. Age Settings


After you've confirmed your height, the system moves on to your age settings. As shown in the figure below, the number flashes on the display. By holding down the button ◀ or ▶ the user sets their age. Press the „SET“ button to confirm the setting.

5. Setting Goals

Once you have confirmed your age, the system will guide you to the goals setting. As shown in the figure below, the number flashes on the display. By holding down the button ◀ or ▶ the user sets a specific target. Press the „SET“ button to confirm the setting.

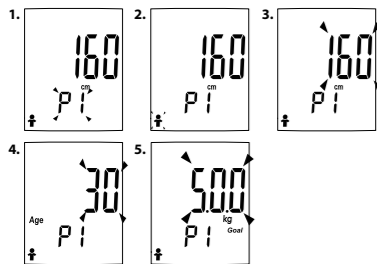LCD display:

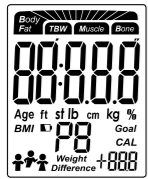BF	Body Fat Analysis Result	%	Percentage
TW	Total Body Water Analysis Result	♂	Male
MM	Muscle Mass Analysis Result	♂	Male Athlete
BM	Bone Mass Analysis Result	♀	Female
Age	Age	♀	Female Athlete
ft	Foot	PB	User ID (P1-P8)
st lb	Stone	CAL	Calorie Analysis Result
lb	Pound	BMI	Body Mass Index
cm	Centimeter	⬇	Low Battery
kg	Kilogram	Goal	Goal Weight
--	Less than the Goal Weight	+	More than the Goal Weight

Package contents:

- Your Truelife FitScale W3 comes with:
1. Truelife FitScale W3
 2. 3 x 1.5V AAA battery
 3. User Manual + Safety Instructions

Basic specifications:

Dimensions: 32 cm x 32 cm x 2.4 cm
 Weight: 1.6 kg
 Display: LCD Digital
 Battery: 3 x 1.5V AAA
 Weighing Range: 5 kg ~ 180 kg
 Margin of Error: 0.1 / 0.2 kg

Manufacturer's Responsibility:

Truelife products are covered by a warranty for manufacturing defects. It is the responsibility of elem6 s.r.o. to repair or replace the defective part or product for such defects or warranty for a comparable part or product at its own discretion. Conversion, alteration or other unauthorized manipulation of the product means exceeding the scope of this warranty. Even if defects, negligence or other errors in the product result in loss or damage. The warranty does not apply to damage resulting from use, accident or normal wear and tear. Elem6 sro assumes no responsibility for any accidents, injuries, death, loss or other requirements relating to or arising out of the use of this product. Under no circumstances is elem6 sro liable for incidental or consequential damages in connection with or resulting from the use of this product or any of its components. Repairment or replacement of spare parts and / or products may be subject to transportation, handling, replacement and / or refilling charges.

SAFETY WHEN WORKING WITH THE DEVICE:

1. Carefully read the enclosed Safety Instructions.
2. For the correct handling of the product, read the most up-to-date version of the manual, available for download at www.truelife.eu.
3. Print errors are reserved.

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