# BabyScale B5 BT

---

**EN - MANUAL**

**DE - ANLEITUNG**

**CS - NÁVOD**

**SK - NÁVOD**

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**SL - NAVODILA**

**HR - UPUTE**

---

## EN – MANUAL

---

### INTRODUCTORY INFORMATION

- Before switching on the scale for the first time, insert the 3 included AAA batteries into the compartment at the bottom of the scale. Make sure they are oriented according to the diagram in the compartment.
- Place the scale on a flat, solid surface or a mat in a dry environment.
- If you use the scale on a soft surface or carpet, the product will not be able to measure the child's weight optimally.
- Do not exceed the maximum weight (20 kg).
- When not in use, do not place anything on top of it or put unnecessary strain on it.
- Do not shake the scale and avoid impacts.
- Avoid contact with water.
- For the rest of the information regarding the safe use of the device, see the Safety Instructions in the box.

### FIRST USE, WEIGHING WITHOUT THE APP

1. Press the **ON/OFF/TARE** button, the scale will turn on, and the display will light up and show „Hi“.
2. When the text disappears and the scale display shows 0, the scale is ready for use.
3. Place the child on the scale, which will then display the measured weight.
4. The measured value on the display flashes 3 times and then remains displayed.
5. Turn off the scale by holding down the **ON/OFF/TARE** button, or the scale will automatically turn off after 2 minutes.

### FIRST USE, WEIGHING WITH THE APP

- The scale can also be used in conjunction with a smartphone app to help you monitor your child's development and growth using measurement history and graphs.

#### 1) Download the app and register

- a. Scan the QR code or enter the URL link and install the app.
- b. Fill in your personal information and register.

[truelife.eu/58605/app](https://truelife.eu/58605/app)

#### 2) Pair the scale with the app

- a. Turn on the app and activate Bluetooth on your smartphone.
- b. In the app, click the **Add device** button or the + in the upper right corner.
- c. From the menu select **Baby scale (Bluetooth)**.

- d. Turn on the scale – the Bluetooth icon will flash in the upper right corner of the display.
- e. In the app, press **Next**, then select the scale.
- f. Now you can name the scale and choose the room where you will use it. Then select **Start Using**.
- g. Choose the units in which you want the app to display weight and height data and confirm your selection.
- h. The scale is successfully paired and ready to use.

### 3) Weighing with the app

- a. Turn on the app and the scale.
- b. Place the child on the scale. The scale immediately displays the child's weight both on its display and in the app.
- c. The measured value on the display flashes 3 times and then remains displayed, and in automatic HOLD mode, it is directly transferred to the application.
- d. You can reset the scale in the application with **Reset**, or if you use manual HOLD, you can lock and transfer the weight to the app with the **Record** button.
- e. At the bottom, you can see the measurement history. When you select any of the recorded measurements, you will see more detailed information including a graphical display.
- f. You can turn off the scale by holding down the **ON/OFF/TARE** button, or the scale will turn itself off after 2 minutes.

## FUNCTIONS

### 1. HOLD function (automatic)

- The TrueLife BabyScale B5 BT scale has a HOLD function, which helps parents accurately weigh restless children who move around on the scale.
- When the baby is placed on the scale, the function is activated and the scale measures the baby's weight when they are moving the least and locks the value on the display or writes it directly into the app.
- Automatic is the default setting of the HOLD function.

### 2. HOLD function (manual)


- If you prefer not to use the HOLD function (or to use it manually), it is possible to switch from the default automatic HOLD to manual HOLD on the scale at any time.
- Switch between automatic and manual HOLD modes by holding down the HOLD button after turning on the scale. When switching to manual HOLD the scale will display OFF, when switching back to automatic HOLD the scale will display ON.
- If you are in manual HOLD mode, the scale will not lock the value on the display. If you wish to lock the scale and transfer it to the application, press the HOLD button on the scale or the **Record** button in the application.

## TARE function

– The TrueLife BabyScale B5 BT has a TARE function, which allows you to weigh your baby's clothes, for example, and then weigh the baby with clothes. To activate the function, follow these steps:

1. Turn on the scale and wait until the display shows 0.
2. Place the clothes on the scale.
3. Press the TARE button – the TARE icon appears on the display and the scale is reset.
4. Remove the clothes from the scale.
5. Place the child on the scale with clothes or wrapped in a blanket and the scale will show you the actual weight of the child.

## CONTROLS
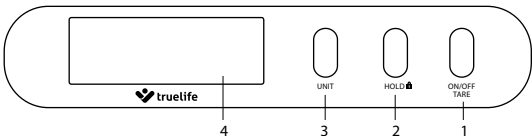
### 1. ON/OFF/TARE

- a. Press to turn on the scale
- b. Press to activate the TARE function
- c. Hold to turn off the scale

### 2. HOLD/🔒

- a. Pressing activates the HOLD function (only in manual HOLD mode)
- b. Holding switches the HOLD function between automatic/manual


### 3. UNIT button

- a. Pressing switches between units of measurement

### 4. Display

## DESCRIPTION OF THE DISPLAYED DATA

1. TARE icon (when the function is turned on)
2. Weight indicator
3. Bluetooth icon
4. Units of measurement
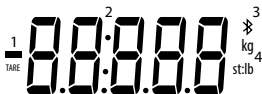
## TROUBLESHOOTING

- If the display shows error message **o-Ld**, the scale is overloaded.
- If the scale displays an error message **Lo**, the batteries need to be replaced.

## BOX CONTENTS

- 1× TrueLife BabyScale B5 BT
- 3× AAA batteries
- 1× Manual
- 1× Safety instructions

## BASIC SPECIFICATIONS

- Dimensions: 67 cm × 33 cm × 5.2 cm
- Weighing range: 50 g ~ 20 kg
- Weighing deviation: 5 g
- Battery: 3× 1.5V AAA
- Operating temperature: 5–35 °C
- Storage temperature: -5–45 °C

## SAFETY WHEN WORKING WITH THE DEVICE

1. Carefully read the enclosed safety instructions.
2. For the correct handling of the product, download and read the latest version of the manual, which can be downloaded at: [www.truelife.eu](http://www.truelife.eu)

**Printing errors and changes to information reserved.**

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