# AquaFloss Compact L100

---

**EN – USER MANUAL**

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**HR – UPUTE ZA UPOTREBU**

---

## EN – USER MANUAL

---


### BOX CONTENTS

- TrueLife AquaFloss Compact L100
- 2× Standard nozzle heads, 1× Dental Plaque nozzle head, 1× Orthodontic nozzle head, 1× Periodontal nozzle head, 1× Tongue Cleaner
- Leatherette case
- USB charging cable
- Manual and safety instructions

### PRODUCT DESCRIPTION

1. Nozzle head
2. Button to release nozzle head
3. On/Off button
4. Display
5. Mode change button
6. Handle
7. Water tank
8. Water tank cap
9. Charging Cable

**Warning:** *There may be residual water in the device. This is distilled water used during quality control testing.*
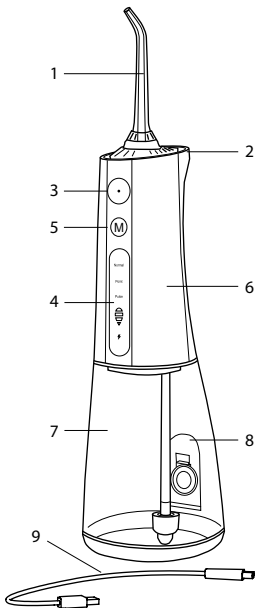
## INSTRUCTIONS FOR USE

Before the first use, remove all packaging material.

If the water flosser has been stored in a low temperature environment during transport, leave it at room temperature for at least two hours before use.

1. Connect the head to the shower. A click indicates correct insertion. You can choose from the following headers:
  - **Standard Head** – For daily plaque removal from interdental spaces and along gum lines
  - **Periodontal Head** – For cleaning periodontal pockets using a rubber tip
  - **Orthodontic Head** – For removing plaque from brackets and wires of braces and other hard-to-reach places
  - **Dental Plaque Head** – For more effective removal of plaque around implants, crowns and bridges
  - **Tongue Cleaner Head** – For effective cleaning of dirt on the tongue
2. Remove the water tank cap and fill the tank with water; the temperature of the water should not exceed 40 °C. The capacity of the water tank is 300 ml. Then replace the cap, ensuring that water does not spill.
3. Hold down the mode change button and select the desired shower mode. Press the same button to set the pressure level (1-4 levels). You can change the mode or pressure even while showering. Modes:
  - **Normal mode** – suitable for daily cleaning. Pressure Range: 82–128 PSI
  - **Point mode** – suitable for effectively removing food residues between the teeth, thanks to a program of three quick and precise sprays under high pressure followed by a short pause. Pressure Range: 41–68 PSI
  - **Pulse mode** – suitable for massaging teeth and gums. Pressure Range: 97–125 PSI
4. Leave your mouth slightly open to allow water to rinse from your mouth into the sink. Press the On / Off button to start the flosser. Direct the water jet at a 90° angle to the teeth and gums. Slowly run the stream of water around your teeth. Direct the stream of water onto the gums and into the space between the teeth. If you have braces, you can direct the water between the joints of the braces.
5. When you have finished flossing, switch off the device using the On / Off button. Open the water tank cap and pour out any remaining water. Switch the device on again with the On / Off button to drain the water from the pump. Press the On / Off button to turn off the device and close the water tank cap.

**Warning:** *To prevent the build-up of bacteria, ensure that no water remains in the device. If you do not intend to use the device for a longer time, make sure that it is thoroughly dried before storing it.*

## CHARGING

The water flosser is not supplied with a fully charged battery. Before using for the first time, fully charge the water flosser.

Plug one end of the USB cable into the charging socket at the rear of the water flosser. Plug the other end of the USB cable into an adapter (adapter not included). Insert the adapter into the plug socket.

It takes approximately 6 hours to fully charge the water flosser.

**Note:** *To ensure maximum battery life, we recommend charging the water flosser once a week.*

## INSERTING AND REPLACE THE NOZZLE HEAD

Insert the nozzle head into the hole at the top of the water flosser. You will hear a click when the nozzle head has been correctly attached.

To remove the nozzle head from the main body, press and hold the release button and then detach the nozzle head from the water flosser.

The nozzle head rotates 360°.

**Warning:** *Only insert and remove the nozzle head when the water flosser is switched off.*

**Note:** *Nozzle heads are items of personal hygiene and are not covered by warranty.*

## MAINTENANCE

Use only water and neutral detergents to clean the device. Never use caustic or abrasive cleaners as they may damage the device. Do not use water whose temperature exceeds 50 °C. Keep the device away from high temperatures and avoid direct sunlight.

### 1. The device

Remove dirt by wiping with a towel or cloth with liquid soap and water, then ensure all water is rinsed off and wipe dry.

Do not immerse the device in water.

**Warning:** *Water must not enter the USB socket.*

### 2. Water tank

Wash the tank with water and wipe with a dry cloth.

Dry the tank thoroughly if you do not intend to use it for a longer time.

**Note:** *For greater ease of cleaning, the water tank can be removed from the body of the water flosser.*

### 3. Nozzle head

Wash the nozzle head with water and wipe with a cloth.

Replace the nozzle head regularly, at least once every three months, or more often depending on wear.

## TROUBLESHOOTING

| Problem                                      | Possible cause  | Solution   |
|--|---|--|
| Water pressure is very low                   | The nozzle head is unclean / battery is low                   | Replace the nozzle head with a new one / recharge the water flosser              |
| No water at all emits from the water flosser | Water tank is empty / The device is tilted at the wrong angle | Refill the water tank / Move the device to the correct position (i.e., vertical) |
| Nozzle head cracks                           | Nozzle head has been used for too long                        | Replace the nozzle head with a new one   |
| Water tank leaks                             | Water tank cap is damaged                                     | Contact the manufacturer to replace the water tank                               |

## SPECIFICATIONS

|                    |              |
|--------------------|--------------|
| Energy consumption | 10 W         |
| Voltage            | 5 V          |
| Battery capacity   | 1500 mAh     |
| Water resistance   | IPx7         |
| Pressure range     | 30 – 130 PSI |
| Reservoir volume   | 300 ml       |
| Weight             | 264 g        |
| Dimensions         | 79×73×214 mm |

## PRODUCT LIABILITY

TrueLife products are guaranteed to be free from manufacturing defects. In the event of such defects, elem6 s.r.o. will repair or replace the defective part or product with a comparable part or product at its sole discretion. The warranty will be deemed void if the device is altered or tampered with. The warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss, or other claims related to or resulting from the use of this product. elem6 s.r.o. shall in no event be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Return or exchange of spare parts and / or products may be subject to shipping, handling, and replacement or replenishment fees.

This product should be stored out of the reach of children. Use by children should take place only under the direct supervision of an adult.

## SAFETY GUIDELINES WHEN OPERATING THE DEVICE

1. Read the enclosed safety instructions carefully.
2. The latest version of the manual, including instructions for operating the device, can be downloaded at: [www.truelife.eu](http://www.truelife.eu).
3. Printing errors and changes of information reserved.

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