
# BELLE IPL E5

---

- EN - USER MANUAL**
- DE - BENUTZERHANDBUCH**
- CZ - UŽIVATELSKÝ MANUÁL**
- SK - POUŽÍVATEĽSKÝ MANUÁL**
- PL - INSTRUKCJA UŻYTKOWNIKA**
- HU - HASZNÁLATI ÚTMUTATÓ**
- SL - UPORABNIŠKI PRIROČNIK**
- HR - UPUTE ZA UPOTREBU**

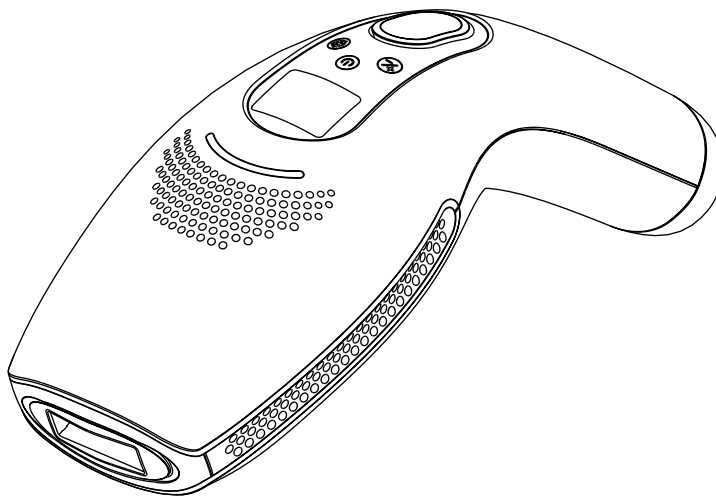
## EN - USER MANUAL

CONTENTS	
BOX CONTENTS .....	1
PRODUCT DESCRIPTION .....	1
SAFETY PRECAUTIONS .....	2
Contraindications .....	2
HOW DOES THE EPILATOR WORK?.....	3
IPL technology.....	3
What can you expect? .....	3
After application .....	3
6th to 12th week of use.....	3
After 12 weeks of use .....	3
HOW TO USE THE EPILATOR? .....	3
Before using for the first time.....	3
Recommended skin colour .....	3
Preparing the skin .....	3
Preparing the device .....	3
EPILATION .....	4
Modes .....	4
Treatment plan .....	4
Lamp head types .....	5
Replacing the lamp head.....	5
AFTERCARE .....	5
DEVICE MAINTENANCE .....	5
CLEANING AND STORAGE .....	5
SIDE EFFECTS .....	5
TROUBLESHOOTING .....	5
TECHNICAL SPECIFICATIONS .....	6
FURTHER INFORMATION .....	6
EU DECLARATION OF CONFORMITY .....	6

## BOX CONTENTS

- TrueLife Belle IPL E5
- HR lamp, AC lamp, SR lamp
- Charging base
- Safety glasses
- Manual

## PRODUCT DESCRIPTION

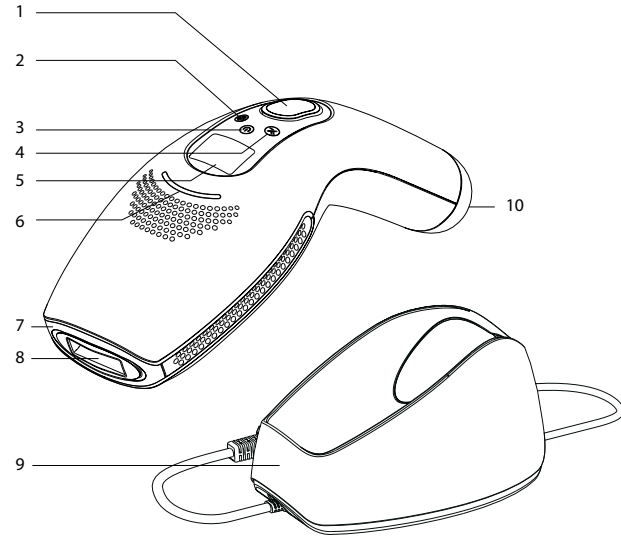1. Treatment button
2. Cooling mode button
3. Power button
4. Level button
5. LCD display

6. Indicators
7. Replaceable lamp head
8. Treatment window
9. Charging base
10. Base cable input

## SAFETY PRECAUTIONS

- To ensure your own safety, do not use the controls on this unit while driving.
- The device must always be used with safety goggles.
- Always keep the device and the adapter dry.
- To reduce the risk of electric shock, do not touch any parts if the unit is broken.
- Do not use the device in a humid environment - for example with running water or a full bath.
- The device may be used by children and persons with reduced physical, sensory or mental capabilities, provided they have been given supervision or instruction concerning use of the device and understand all risks associated with the use of the device.
- The device must not be used by children under 15 years of age. Persons aged 15-18 only under the supervision of a guardian.
- The device is equipped with a transformer and can only be used with the original one supplied.
- Always check the device before use. Do not use the device if any part is damaged.
- If the adapter is damaged, do not use it. There is a risk of electric shock.
- The light emitting treatment lamp/window may be very hot after use. Do not touch the inside of the window immediately after application.
- Before using the device, you should clean the skin and ensure that it is completely hairless, dry and free of greasy substances.
- Make sure that nothing obstructs the flow of air from the vents of the device.
- Never subject the device to shocks, shake it or drop it.
- If you move the device from a very cold environment to a very warm environment or vice versa, wait about 3 to 4 hours before using it.
- Store the device in a dry place where it will not get covered in dust.
- Do not expose the device to direct sunlight or UV radiation, otherwise it may be damaged.
- Never leave the device unattended when it is switched on. Always switch off the device after use.
- For hygienic reasons, the device is intended for one user.
- The device is only intended for removing unwanted hair from parts of the body from the face and neck down. Do not use it for any other purpose.
- The device is not washable. Never immerse the device in water or wash it under the tap.
- The device is not intended for children. It is not a toy and children should not use it.
- Do not mark the treated areas with a pencil or pen. It could cause skin burns.
- Use only with settings that are appropriate for your skin type. When used with settings other than those recommended, the risk of skin reactions and side effects may increase.
- Never leave the device unattended when it is switched on. Always switch off the device after use.
- Do not use the device if any of the cases listed in the contraindications section apply to you.
- Do not treat the same area of skin more than once during one session. This does not increase the effectiveness of the treatment, but increases the risk of skin reactions.
- Hair removal using lasers or intense pulsed light sources can cause increased hair growth in some individuals. According to available data, the highest risk of this reaction is in women of Mediterranean, Middle Eastern and South Asian descent treated on the face and neck.
- The box contains small pieces that can be dangerous for children. Always keep the product

out of the reach of children. The bags can cause suffocation if placed on the head and the pieces they contain are a choking hazard.

- For pluggable equipment, the socket-outlet shall be installed near the equipment and shall be easily accessible.

## Contraindications

- Never use the device if you have skin type VI (for type V at your own discretion). In this case, you are at high risk of skin reactions, such as excessive or insufficient pigmentation, significant redness or burns. For a list of skin types, see page 3.
- Never use the device if you are pregnant or breastfeeding. The device has not been tested on pregnant or lactating women.
- Never use the device if you have any active implant, such as a pacemaker, neurostimulator, insulin pump, etc.
- Do not use the device if you are taking any of the following medicines:
  - if your skin has been treated in the last week with alpha hydroxy acids (AHA), beta hydroxy acids (BHA), topical isotretinoin and azelaic acid.
  - if you have taken any form of isotretinoin Accutane or Roaccutane in the last six months. This treatment thins the skin so that it can then be more prone to cracks, injuries and irritation.
  - If you are taking photosensitive products or medicines, check the package leaflet of the medicine and never use the device if the medicine causes photoallergic reactions, phototoxic reactions or if you should avoid the sun while taking this medicine.
  - if you are taking anticoagulants (including large doses of aspirin).
  - if you are taking analgesics that reduce your skin's sensitivity to heat.
  - if you are taking immunosuppressive medicines.
- Never use the device:
  - if you have received radiation therapy or chemotherapy in the last 3 months.
  - if you have undergone surgery in the last three weeks in the areas you want to treat.
- Do not use the device if you have or have had the following illnesses:
  - if you have diabetes or other systemic or metabolic diseases.
  - if you have had heart failure.
  - if you have a light-sensitive disease such as polymorphic light flare (PMLE), solar urticaria, porphyria, etc.
  - if you have had a collagen disorder, including keloid scarring or poor wound healing.
  - if you suffer from light-sensitive epilepsy.
  - if you have light-sensitive skin that easily get a rash or have an allergic reaction.
  - if you have a skin condition, such as active skin cancer, a history of skin cancer, or any other localized outbreak of cancer in the areas you want to treat.
  - if you have a history of a venous disorder, such as varicose veins or spider veins in the areas you want to treat.
  - if you have any bleeding disorder.
  - if you have a history of immunosuppressive disease (including HIV infection or AIDS).
- Do not use the device:
  - if you have inflammation, eczema, burns, inflamed hair follicles, open wounds, abrasions, herpes simplex virus, injuries or lesions and hematomas in areas you want to treat.
  - If you have irritation (redness or cuts), sunburn, recent sunburn or artificially tanned skin.

- Do not use the device on the following places/areas:
  - birthmarks, freckles, dilated veins, darker pigmented areas, scars, skin irregularities without consulting a doctor. Application to these areas can lead to burns and skin discoloration, which could potentially make it more difficult to identify skin diseases.
  - warts, tattoos or permanent makeup.
  - around or near the eyes and eyebrows.
  - lips, nipples, areolas, small labia, vagina, anus and inside the nose and ears.
  - men must not use the device on the face and neck, including all areas where beards grow, or the entire genital area.
  - in areas where you use long-acting deodorants. Skin reactions may occur.
  - on or near artificial objects, such as silicone implants or piercings.

## HOW DOES THE EPILATOR WORK?

### IPL technology

Epilation is performed using intense pulsed light. The whole process is painless, leaving your skin smooth for up to several months.

Using IPL technology, the epilator emits gentle pulses of light in the wavelength range of 510 to 1100 nm below the surface of the skin. This targets the melanin in the hair follicle and helps break the cycle of hair re-growth. Light pulses create heat that warms the hair follicle, destroying it which leads to the hair falling out.

### What can you expect?

#### After application

You should not see any significant side effects of the treatment immediately after application. During the first weeks of use, hair growth will still be visible. These are probably hairs that were not in their growth phase during application, when IPL technology is most effective.

#### 6th to 12th week of use

You should already see a reduction in hair growth. However, many hairs will still not have been treated in their growth phase. It is important to continue weekly treatment.

#### After 12 weeks of use

Approximately after completing the 12-week plan, you should already notice a significant reduction in hair in the treated area. The remaining hairs should be finer and lighter. Continued monthly treatment or treatment as needed should keep your skin hairfree.

## HOW TO USE THE EPILATOR?

### Before using for the first time

#### Recommended skin colour

TrueLife IPL epilators work effectively and safely on different skin tones. In general, IPL epilators are not recommended for dark skin and light hair.

To determine if an IPL epilator is right for you, use the Fitzpatrick scale, which distinguishes 6 degrees of skin colour, and determine if an IPL epilator is right for you.

Type I		Pale skin, freckles, blond or red hair, blue eyes. Very sensitive to the sun. It always burns, never tans.	Suitable
Type II		Light skin, blond or red hair, blue, green or brown eyes. Very sensitive to the sun. It usually burns, tans minimally.	Suitable
Type III		Creamy fair skin, any hair and eye colour, common type. Sun sensitive skin. Sometimes it burns slightly, tans evenly.	Suitable
Type IV		Slightly brown skin, typically the Mediterranean olive skin tone. It rarely burns, tans well.	Suitable
Type V		Dark brown skin, typically Middle Eastern or Latin American skin tone. It rarely burns, tans very easily.	Suitable but with caution
Type VI		Deeply pigmented dark brown to black skin. It never burns, it tans very easily.	Unsuitable

### Preparing the skin

Before use, make sure that your skin is in good condition without injury. It must also be shaved, clean and dry (without any residue such as creams, make-up, antiperspirants or deodorants).

Before the first epilation, we recommend performing a patch test on a small area of skin. Perform the test 24 hours before epilation to check for any side effects.

#### How to test

1. Choose an area of skin close to the area you want to epilate.
2. Switch on the device and select level 1.
3. Place the device on the skin, press the treatment button and the device will emit a flash. Emit two flashes like this, corresponding to two treatment windows on the surface.
4. Increase the setting by one level and repeat the flashes. Repeat for the other levels until the flashes are uncomfortable. Once you reach a level that is uncomfortable for you, end the test. Use the last level, which is pleasant for you, during the subsequent epilation.
5. After the test, wait 24 hours and check for a skin reaction. If the skin looks normal, you can proceed with the first epilation. A patch test should be performed on all areas you are going to treat.

### Preparing the device

1. Remove the device and its accessories from the box.
2. Before starting treatment, make sure that the device is clean, dry and undamaged. You can clean the epilation window.
3. Connect the charging base to an electrical outlet and then the device to the charging base.


4. The device is now ready for use.

*Note: Before switching on, make sure that the lamp head is installed correctly. The epilator cannot be switched on without the lamp inserted.*

## EPILATION

If the patch test is successful and your skin is prepared, you can start epilating.

1. Press the power button for 3 to 5 seconds to turn on the device.
2. The LED indicator starts flashing blue and white. When the device is ready for use, the indicator lights up white only.

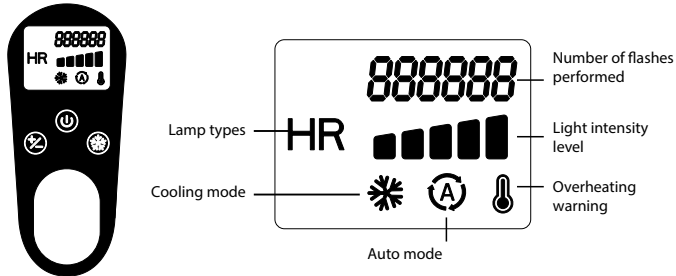
3. On the display you can see: the type of lamp, the number of flashes performed, the light intensity level, the cooling mode icon, the auto mode icon and the overheating warning icon.
4. Press the level button to set your desired light level (level 1 lowest, level 5 highest), which you determined during the patch test.

*Note: When the device is switched on, the light intensity level is set to level 1 by default.*

5. When you have everything ready, touch the skin where you want to epilate. As shown in the picture, the indicator lights up blue and the device is ready to emit light.

*Note: The device must be perpendicular to the epilated area.*

6. Press the treatment button to emit a flash. The flash should make you feel warm.
7. Place the device on the next area you want to epilate. After each flash, it takes up to 3 seconds for the device to be ready for the next flash. You can emit a flash when the blue indicator lights up again.
8. To ensure complete treatment of all areas, emit flashes close together. Effective light comes only from the light emitting treatment window, so the treated areas need to overlap. However, always be sure to send the flash to the same area only once. Irradiation of the same area does not double the effectiveness of the treatment, but increases the risk of skin reactions.
9. When you have finished epilating, hold down the power button for 3 to 5 seconds and the device will switch off.

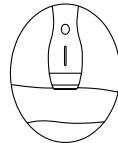10. Disconnect the device from the power supply.

*Note: The device may enter overheat protection mode after a long period of operation. The device cannot be used in this mode until its temperature drops. The device can be used again as soon as the overheating warning icon disappears.*

*Note: If redness or burning occurs for several minutes, soothe the skin by applying a coolant or a cotton pad.*

*Note: If you notice any discoloration on the skin (such as fading or darkening), stop using the device immediately and undergo a medical examination to determine the cause.*

## Modes

The device offers two epilation modes for convenience. In addition to the standard single flash mode, it is also possible to use automatic mode.

### Standard mode

Standard mode is suitable for treating small and hard-to-reach areas. These include the face, armpits, bikini area, knees, shins and ankles.

### Auto mode

Auto mode offers convenient and fast use on larger areas such as feet.

Hold down the level button for 3 seconds. The device switches to automatic mode. Move the device slowly and smoothly from the bottom to the top of the epilated area. Make sure that the device is in constant contact with the skin so that it can emit flashes at all times.

### Cooling mode

Cooling mode is automatically switched on. During epilation, the skin is cooled at the same time making the treatment less painful.

If you want to turn off the automatic cooling when epilating, press the cooling mode button. Pressing the button again to switch it on.

Cooling mode can also be used on its own. Press the cooling mode button for 3-5 seconds. You can then place the device on the area and slowly slide it to cool the area. To turn off cooling mode, press the cooling mode button again for 3-5 seconds.

## Treatment plan

The hair growth cycle varies from person to person. It ranges around 16 months depending on which part of the body. Due to the fact that only hairs in the growth phase respond to the epilation, several treatments are needed to achieve permanently smooth skin. For these reasons, we recommend the following plan.

### Initial phase

For the first 4 to 12 weeks, we recommend epilating once a week.

If hair growth stops before, it is not necessary to continue treatment until the end of week 12. In this case, you can move on to the next stage.

### Maintenance phase

After the initial phase, we recommend treatments every 4 to 8 weeks when you notice hair growth. This phase helps to preserve the results achieved and maintain smooth skin for many months. The time between treatments may vary with regard to individual hair growth and also to each specific area of the body.

## Lamp head types

The device can be used with three different lamp heads. Each lamp has a different function.


**HR lamp** – The HR lamp is designed for hair removal.

**SR lamp** – The SR lamp is designed to rejuvenate the skin.

**AC lamp** – The AC lamp is designed to remove acne.

## Replacing the lamp head

The device allows you to remove and replace the lamp head. Follow the instructions below.
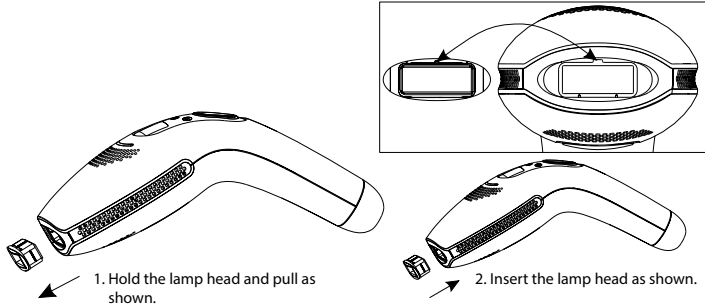
### How to use an AC lamp for acne

1. Before each use of the epilator, first squeeze out any blisters and sores on your face and cleanse your skin.
2. Wear the protective goggles provided in the package.
3. Turn on the epilator and make sure it has the AC lamp on.
4. Apply the epilator with the AC lamp on to the acne affected skin and flash the button. Apply the flashes to the individual acne areas. Give a maximum of 3-5 flashes per spot. Always take short breaks between flashes.
5. Repeat this process every 2-3 days.

After the 10th treatment, you can stop using the facial epilator - your skin should be visibly healthier and without signs of acne damage.

### How to use the SR lamp for skin rejuvenation

1. Before each use of the epilator, cleanse your skin first.
2. Put on the protective goggles provided in the package.
3. Turn on the epilator and make sure it has the SR lamp on.
4. Place the epilator with the SR lamp on your skin and use the button to release the flash. Apply the flashes to individual areas of dry skin. Give a maximum of 3-5 flashes per spot. Always take short breaks between flashes.
5. Repeat this process once every 3 days for dry skin, and once a week for normal skin.

## AFTERCARE

- After epilating, you can safely apply lotion, cream, deodorant, moisturizer or other cosmetic products to the treated area.
- If you experience skin irritation or redness after treatment, wait until these symptoms have disappeared before applying any product to the skin. If you notice skin irritation or redness after applying the product, rinse it with water.
- After epilating, we recommend not exposing the epilated area immediately to sunlight, or to protect it with a product with a protective factor. It is also not advisable to visit tanning booths and solariums within 24 hours after epilation.

## DEVICE MAINTENANCE

- Maintaining this device ensures trouble-free operation and reduces the risk of damage.
- Keep the device away from extreme humidity and temperatures and do not expose the device to direct sunlight or ultraviolet light for extended periods of time.
- Never use compressed air, scouring pads, abrasive or aggressive cleaning agents such as petrol or acetone to clean the device.
- Do not expose the device to sudden and extreme changes in temperature. This can cause moisture to condense inside the device, which could damage the device. In the event of moisture condensation, allow the device to dry completely before using it again.
- Never attempt to disassemble, repair or modify the device in any way. Disassembly, modification or any attempt at repair will void the warranty and may result in damage to the device or even injury or property damage.
- Do not store or carry flammable liquids, gases, or explosive materials in the same place as the device, its parts, or accessories.
- Overheating can damage the device.

## CLEANING AND STORAGE

- After epilation, switch off the device and let it cool down. Be sure to unplug it. Inspect all parts for damage and then wipe them with a cloth or other lint-free material - it can be slightly damp. Do not wash or immerse the device in water or other liquids.
- Store the device in a dry place where it will not collect dust, at a temperature from 0°C to 40°C. Protect the treatment window and skin sensor from damage.

## SIDE EFFECTS

- Depending on your skin type, using the device may be slightly uncomfortable. It may cause some heat or reddening of the skin. This is a normal.
- Adverse reactions usually occur immediately after application or within 24 hours. In rare cases, they may appear later. If you get an unusual reaction, talk to your doctor.

## TROUBLESHOOTING

If a malfunction occurs while using the device, follow the instructions below to find a solution. If you do not find the answer to your problem here, please contact the manufacturer's customer service.

Problem	Solution
Cannot start up.	Make sure the adapter is connected properly. Check that the lamp is inserted in the epilator. The epilator will not turn on without the lamp inserted.
Flashes do not work.	Make sure that the light-emitting treatment window is in contact with the skin properly and that the lights are blue.
Constant blue indicator light.	Make sure that the lamp is not located in an area where there are metal parts or the effects of inductive capacity.
Device does not work in small areas.	Use the correct lamp head for small areas or areas that are difficult to access. It is necessary to ensure that the skin sensors are in contact with the skin.
Device is broken, cracked or does not visually look right.	If you have any doubts about the safe use of the device or damage to the device, do not use it and contact the manufacturer.
The LCD display shows E0.	Something happened to the fan. After the device cools down, turn it off and on again.
The LCD display shows E1.	Voltage problem. Contact the manufacturer.
The LCD display shows E2.	Improperly installed lamp head. Remove the lamp and replace it according to the manual.

## TECHNICAL SPECIFICATIONS

Input voltage	100-240V
Input frequency	50/60 Hz
Output voltage	24V
Output current	4A
Output power	96W
Average efficiency in active mode	88.55%
Efficiency at low load (10%)	91.93%
Power consumption in the no-load state	< 0.15W
Frequency	50/60 Hz
Operating temperature	5 - 35°C
Storage temperature	0 - 40°C
Storage humidity	25% - 75% RH
Operating pressure	860 - 1060 Pa
Emitted wavelengths	510 - 1100 nm
Optical light	2.8 - 4.8 J/cm <sup>2</sup> according to the selected setting
Maximum optical energy	14.88 J
Pulse duration	1 - 3 ms
Number of levels	1 - 5
Number of flashes	600 000

## FURTHER INFORMATION

1) For households: The symbol (crossed-out wheelee bin) on the product or in the accompanying documentation means that used electrical and electronic products must not be disposed of along with your other household waste. In order to dispose of the product properly, please return it to a designated collection point, where it will be accepted free of charge. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. Contact your local authority or the nearest collection point for further details. Improper disposal of this type of waste may result in penalties in accordance with national regulations. Information for users on disposal of electrical and electronic equipment (corporate and business use): For proper disposal of electrical and electronic equipment, please contact your dealer or supplier for detailed information. Information for users on disposal of electrical & electronic equipment in other countries outside the European Union: The above symbol (crossed-out wheelee bin) is only valid in European Union countries. For proper disposal of electrical and electronic equipment, request detailed information from your authorities or equipment dealer. Everything is expressed by the crossed-out wheelee bin symbol on the product, packaging, or printed materials.

2) Please contact your dealer for warranty repairs. In case of technical problems and questions, please contact your dealer who will inform you of the next steps to be taken. Observe the rules for working with electrical equipment. The user is not authorised to disassemble the equipment or replace any of its components. There is a risk of electric shock if the covers are opened or removed. Incorrect assembly and reconnection of the equipment also exposes you to the risk of electric shock.

The warranty period for the products is 24 months unless otherwise stated. The warranty does not cover damage caused by non-standard use, mechanical damage, exposure to harsh conditions, handling contrary to the manual and normal wear and tear. The warranty period is 24 months for the battery and 6 months for its capacity. For more information on the warranty, please visit [www.elem6.com/warranty](http://www.elem6.com/warranty)

Neither the manufacturer, importer nor distributor shall be liable for any damage caused by the installation or improper use of the product.

## EU DECLARATION OF CONFORMITY

elem6 s.r.o. hereby declares that the TrueLife IPL BELLE E5 complies with the essential requirements and other relevant provisions of Directive 2014/30/EU and 2014/35/EU. All TrueLife products are intended for sale without restriction in Germany, the Czech Republic, Slovakia, Poland, Hungary and other EU member states. The full Declaration of Conformity can be downloaded from [www.truelife.eu/support/doc](http://www.truelife.eu/support/doc)

## Manufacturer:

elem6 s.r.o.,  
Braškovská 15, 16100 Prague 6, Czech Republic

Printing errors and changes in manual and safety instructions reserved.

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